

Extending PVSCAN to Meet the Market Needs for High-Speed, Large-Area Scanning

B. Sopori, W. Chen, Y. Zhang, T. Hemschoot,
and J. Madjdpour

*Presented at the 9th Workshop on Crystalline Silicon
Solar Cell Materials and Processes
Breckenridge, Colorado
August 9-11, 1999*



NREL

National Renewable Energy Laboratory

1617 Cole Boulevard
Golden, Colorado 80401-3393

NREL is a U.S. Department of Energy Laboratory
Operated by Midwest Research Institute • Battelle • Bechtel

Contract No. DE-AC36-98-GO10337

NOTICE

The submitted manuscript has been offered by an employee of the Midwest Research Institute (MRI), a contractor of the US Government under Contract No. DE-AC36-99GO10337. Accordingly, the US Government and MRI retain a nonexclusive royalty-free license to publish or reproduce the published form of this contribution, or allow others to do so, for US Government purposes.

This report was prepared as an account of work sponsored by an agency of the United States government. Neither the United States government nor any agency thereof, nor any of their employees, makes any warranty, express or implied, or assumes any legal liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. Reference herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States government or any agency thereof. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States government or any agency thereof.

Available electronically at <http://www.doe.gov/bridge>

Available for a processing fee to U.S. Department of Energy
and its contractors, in paper, from:

U.S. Department of Energy
Office of Scientific and Technical Information
P.O. Box 62
Oak Ridge, TN 37831-0062
phone: 865.576.8401
fax: 865.576.5728
email: reports@adonis.osti.gov

Available for sale to the public, in paper, from:

U.S. Department of Commerce
National Technical Information Service
5285 Port Royal Road
Springfield, VA 22161
phone: 800.553.6847
fax: 703.605.6900
email: orders@ntis.fedworld.gov
online ordering: <http://www.ntis.gov/ordering.htm>



Printed on paper containing at least 50% wastepaper, including 20% postconsumer waste

Extending PVSCAN to Meet the Market Needs for High-Speed, Large-Area Scanning

Bhushan Sopori, Wei Chen, Yi Zhang, Tess Hemschoot, and Jamal Madjdpour

National RenewableEnergy laboratory
1617 Cole Boulevard
Golden, CO 80401

INTRODUCTION

PVSCAN is a versatile instrument that has many applications in the PV industry, including high-speed mapping of material and cell parameters such as defect density, reflectance, and LBIC response. Recently, the PV community has been interested in acquiring this instrument for material and cell analyses and for process monitoring. NREL has made a commitment to devote internal resources to either build PVSCAN or provide technical assistance for those who wish to license the instrument. Consequently, it has become necessary for us to document the pertinent information about the system and to prepare a suitable operational manual. We have also learnt that each buyer may have somewhat different needs that must be accommodated while making each instrument. In particular, there is increasing interest in upgrading the system to have the capability to scan larger substrates and cells, with higher speeds and finer resolution.

This paper explores various issues that arise in developing a commercial instrument such as PVSCAN. Emphasis is on the technical details of the ability to scan fast and the detrimental effects this fast scan can have on the image quality of various material/cell parameters.

A BRIEF DESCRIPTION OF PVSCAN

PVSCAN uses the optical scattering from a defect-etched sample to statistically count the density of defects. It shines a laser beam on the surface of the defect-etched wafer and measures the (integrated) intensity of the reflected (scattered) light. The total integrated reflected light is proportional to the number of scattering centers [1]. Thus, this system provides a signal that is proportional to the local dislocation density. By scanning over the sample, the instrument can map the defect distribution. Furthermore, grain boundaries and dislocations have different characteristic scattering-patterns, and PVSCAN uses these patterns to distinguish between different kinds of defects, as illustrated in Figure 1. Figure 2 is a schematic of the optical system.

PVSCAN also provides a quantitative means of measuring the LBIC response of a solar cell at two different wavelengths (0.633 and 0.905 μm) of light excitation. This enables the instrument to separate the near-surface and the bulk recombination characteristics of the cell. The photo-current response for each excitation can be measured and saved by the computer as the external

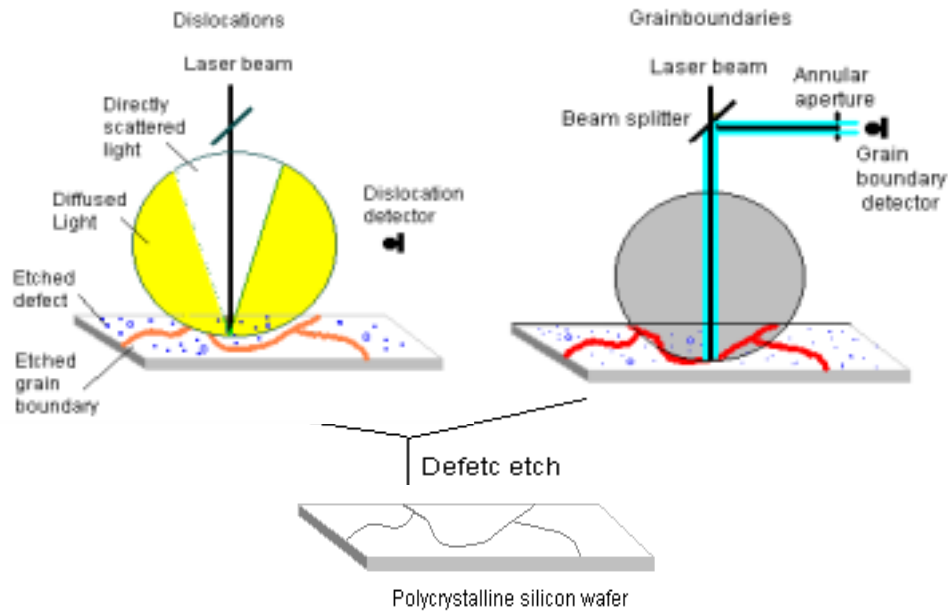


Figure 1. PVSCAN uses differentially reflected light to distinguish between various defect types. After treatment with a special etch, dislocations scatter the incoming laser beam in a cone-shaped pattern which the PVSCAN captures with an integrating sphere (left). Grain boundaries reflect the beam nearly directly, and the PVSCAN redirects a portion of that reflection to a separate detector.

response. PVSCAN's capability of mapping reflectance provides an important step toward identifying losses in the cell because reflectance is a major cause of "efficiency loss" for solar cells. By combining the LBIC reading with reflectance losses, the instrument calculates the cell performance as a function of the light that is actually transmitted into the cell to get an internal photoresponse; this is the core information needed to improve cell performance.

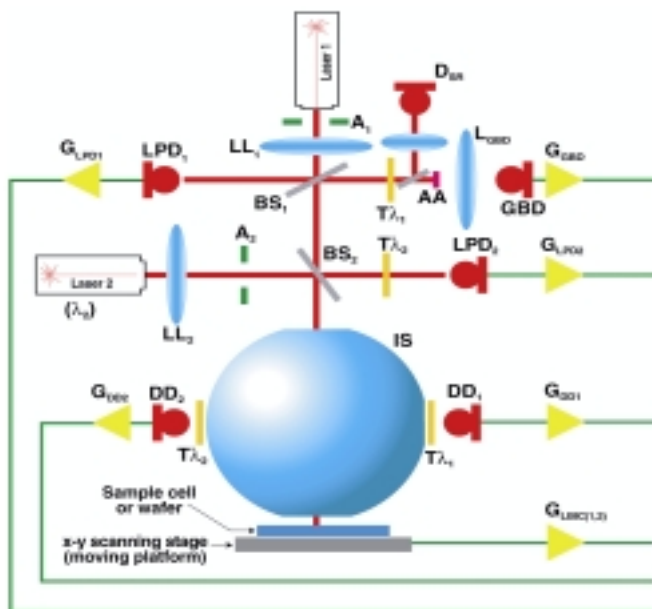


Figure 2. A schematic of the optical system of PVSCAN.

DOCUMENTATION

To build a system, one needs a host of documentation of sufficient detail. The documentation includes:


- a complete set of drawings
- listing all parts/naming the parts
- listing functional modules
- an instruction manual.

NREL is in the process of documenting all details that will allow an outside job-shop to fabricate various parts and the chassis.


MEETING THE TECHNICAL CHALLENGES

1. Larger scan area

There is increasing interest in going to larger wafers and solar cells. The PVSCAN was originally designed to accommodate 4-in. x 4-in. wafers and cells. However, industry is already making commercial cells that exceed this size. We have received requests to increase the ability of PVSCAN to scan 8-in. x 8-in. cells. It is a straightforward matter to achieve this by switching to an X-Y stage with a larger travel, but an increase in the scan area generally leads to an increase in the scanning time (if the resolution is kept the same). In general, an increase in the sample area by a factor of 4 requires an increase in the scan time by a same factor.

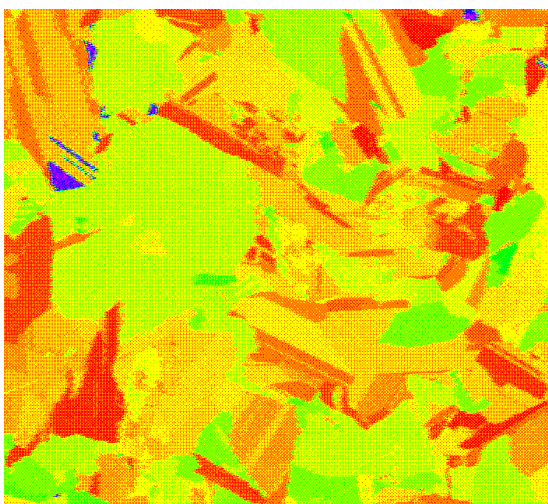

(a) A defect map of a 2-in.x2-in. area at the scanning speed of 3 inch/s. The resolution is 0.001 inch. This scan took about 1 hour.

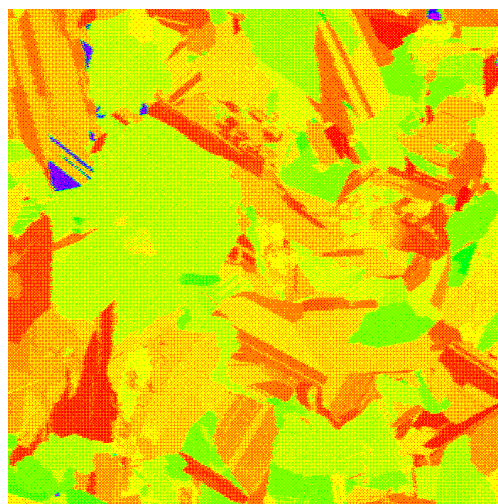

(b) A defect map of the same area at the scanning speed of 1 inch/s. The resolution is 0.001 inch. This scan took about 2.5 hours.

Figure 3. Images of a 2-in. x 2-in. wafer taken at two different scanning speeds, showing nearly identical image quality.

Our approach to increasing the scan area consists of switching to an X-Y stage that has a capability of scanning at a higher speed, typically up to 4 inches per second. Thus, we can scan a 8-in. by 8-in. area, keeping the scanning time and resolution the same as the original PVSCAN. It is crucial that as the speed is increased, the image quality should not deteriorate. Figure 3 compares the quality of the reflectance images of a sample scanned at two different speeds in our system. The resolution is 25 μm , and the difference in the image quality is essentially indistinguishable.

2. Higher resolution

The resolution of a scanning system is a complicated parameter that is determined by:

- Step size of the motors
- Optical beam size
- Time interval between the data-taking events
- Time constant of the electronic signal channel.

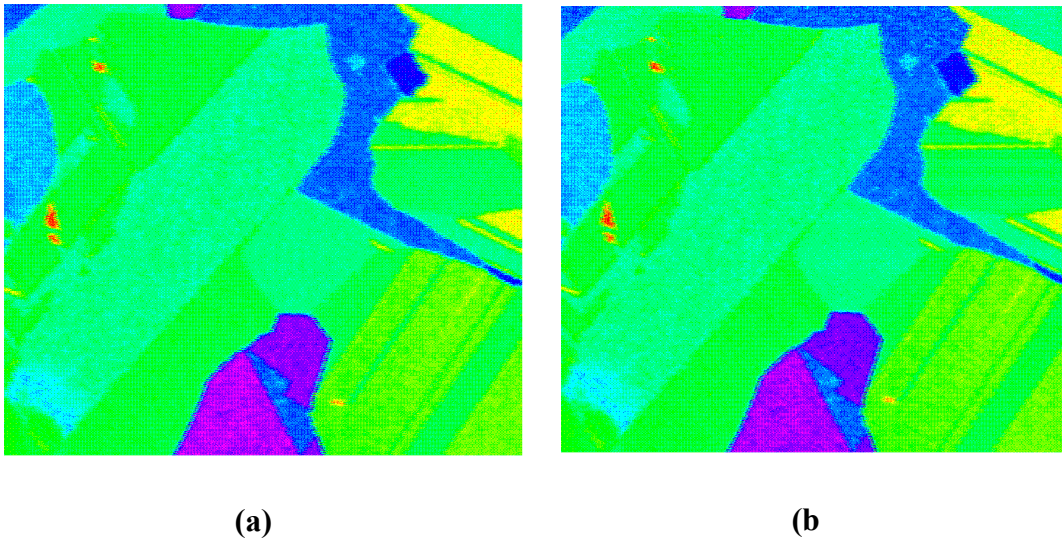


Figure 4. The defect maps of a 0.5-in.x0.5-in. area. The scanning speed is 0.75 inch/s. The resolutions are (a) 0.001 inch and (b) 0.00025 inch.

In general, the PVSCAN uses a beam size of about 200 μm in the defect mapping mode and smaller size in the LBIC mode. For a motor step of 25 μm , a 50 μm spatial resolution in the image can be achieved. Figure 4 shows that with a typical beam size the scanning resolution has little or no effect on the image quality. Even under higher magnification, Figure 5, the effect is primarily observable as a slight loss in the Y-resolution.

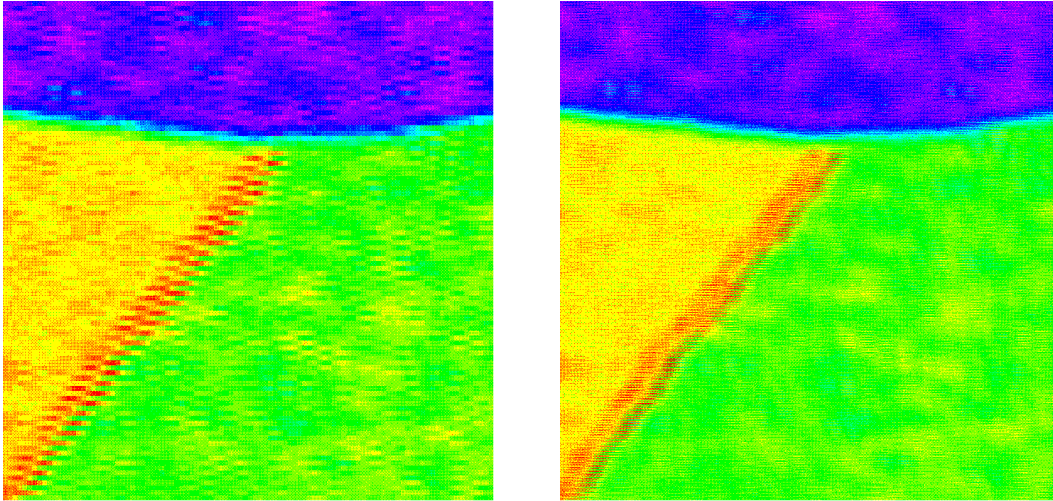


Figure 5. A zoomed-in view of a small region of the sample in Figure 4, showing no effect on the image quality, even when the display is zoomed into a small area.

3. *Effect of beam size*

It is desirable to have beam size about the same as the lateral step size of the motor. But, in practice, a small beam size makes stringent demands on the optical system. Furthermore, it is important that the beam size exceeds a certain value so that statistical counting of defects is valid. In general, the beam-size requirements for most applications of PVSCAN are only moderate. However, it must be pointed out that mathematically the beam and the motor step have the effect of generating a signal that is an overlap integral of the two. This means that the signal in the “forward” and the “reverse” directions will be shifted laterally. An example of this is illustrated in the scans of Figure 6.

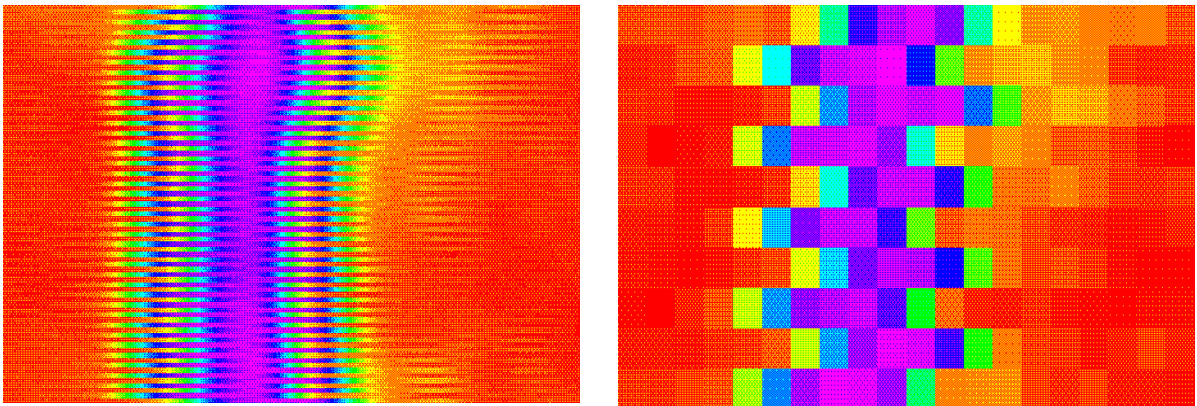


Figure 6. Illustration of the beam-size effect on the scanned image. (a) LBIC image in the vicinity of a metal finger on a solar cell, showing displacement of the signal in the “forward” and the “reverse” directions. (b) A zoomed-in view.

4. Higher-speed scanning

It is generally desirable to have as high a scan speed as possible. However, scanning speed is limited by :

- Maximum speed of the motor(s), resolution of the stepping motor, and the mode of operating the motor (open Vs closed-loop configuration)
- Data-acquisition rate and data-transfer capability
- “Blank” time, in which a Y-step is taken (this time is also used for transferring data from the internal memory to the disk)
- RC time constant of the signal channel(s) and its changes with the gain.

These issues can be optimized to obtain the highest scanning speed. Our system can be operated up to 4 ips without degrading the image quality (as seen in Figure 3).

5. Improvements in the software

The software has been improved to include:

- Statistical module that displays distributions (see Figure 7, showing distribution of defect density in a multicrystalline Si wafer) and other statistical functions
- Additional options on “save” images/scans
- Ease of calibration

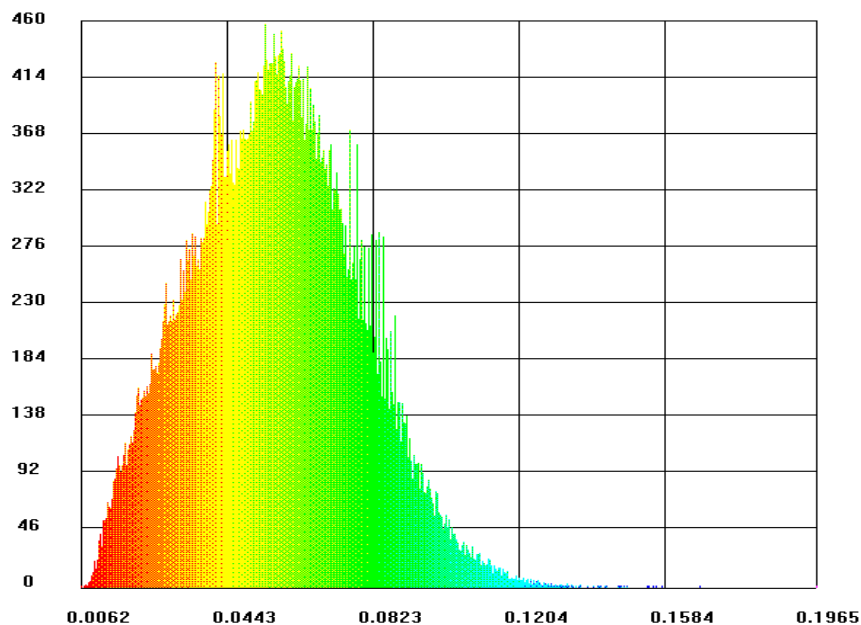


Figure 7. Distribution of dislocations in a mc-Si wafer. The Y-axis shows number of counts and x-axis is the dislocations density, $\text{cm}^{-2} \times 10^4$

CONCLUSION

We have documented the PVSCAN so that it can be built at a job-shop. The operation of the system has been critically analyzed to help the user make an informed decision on setting up the scan parameters and interpret the results. A detailed instruction manual is being written that incorporates examples and results of actual scans. The software has been upgraded to include a number of useful features.

REPORT DOCUMENTATION PAGE			Form Approved OMB NO. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503.				
1. AGENCY USE ONLY (Leave blank)		2. REPORT DATE October 1999		3. REPORT TYPE AND DATES COVERED Conference paper
4. TITLE AND SUBTITLE Extending PVSCAN to Meet the Market Needs for High-Speed, Large-Area Scanning			5. FUNDING NUMBERS	
6. AUTHOR(S) B. Sopori, W. Chen, Y. Zhang, T. Hemschoot, J. Madjdpour			C TA: PV902503	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES)			8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES) National Renewable Energy Laboratory 1617 Cole Blvd. Golden, CO 80401-3393			10. SPONSORING/MONITORING AGENCY REPORT NUMBER CP-520-26932	
11. SUPPLEMENTARY NOTES				
12a. DISTRIBUTION/AVAILABILITY STATEMENT National Technical Information Service U.S. Department of Commerce 5285 Port Royal Road Springfield, VA 22161			12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words) PVSCAN is a versatile instrument that has many applications in the PV industry, including high-speed mapping of material and cell parameters such as defect density, reflectance, and LBIC response. Recently, the PV community has been interested in acquiring this instrument for material and cell analyses and for process monitoring. This paper explores various issues that arise in developing a commercial instrument such as PVSCAN. Emphasis is on the technical details of the ability to scan fast and the detrimental effects this fast scan can have on the image quality of various material/cell parameters.				
14. SUBJECT TERMS			15. NUMBER OF PAGES	
			16. PRICE CODE	
17. SECURITY CLASSIFICATION OF REPORT Unclassified	18. SECURITY CLASSIFICATION OF THIS PAGE Unclassified	19. SECURITY CLASSIFICATION OF ABSTRACT Unclassified	20. LIMITATION OF ABSTRACT UL	